#### Yichao

For shutting off, go through the checklist before turning off the computers.

- Turn off equipments
- Turn off computers

Wait before everyone else finish their equipment shutoff before turning off the computers. For powering up, first turn the computers on before turning on the equipments.

- Turn on computers
- Turn on equipments
- Check computer connection to devices

Also wait before everyone else finish their corresponding steps before starting the next one.

#### Record settings before turning things off.

All the lasers should be turned off.

All the home built  $\pm 15$  V, 5V power supplies should be on UPS. 24 V ones are mainly powering amplifiers and can be shut down.

All locks should be disabled.

All temperature controllers should be on UPS and left on.

For each equipment, check if it is in used. Then check the power and make sure if is/should be on UPS. Make sure the "On UPS" label is correct.

## Equipments are mostly located around the laser table facing the wall and the computer.

### Equipment list

ID	Equipment	Location		
	Laser table facing computer or the wall top -	$\rightarrow$ down, west side $\rightarrow$ south side		
1	Home power $(\pm 15,5)$	Lv2 above 623 nm box		
2	Benchtop power 2+1 Chn	Lv2 above 1.0 Cs MOT facing computer		
3	Benchtop power 1 Chn	Lv2 above 1.0 Cs MOT facing computer		
4	Benchtop power 1 Chn (x3)	Lv2 above 1.0 Cs MOT facing wall		
5	Home power (24) (x2)	Lv2 above 1.0 Cs Raman		
6	Benchtop power for 1.0 Cs Raman	Lv2 above New Focus lock box		

7	High voltage power	Lv2 above timebase drivers		
8	Benchtop power for 1.0 Na Raman seed	Lv2 above oscilloscope		
9	Benchtop power 2+1 Chn for Cs MOT beat lock	Lv1 above 1.0 Cs MOT facing computer		
10	Till driver v2 for 1.0 Cs MOT	Lv1 above benchtop powere supply		
11	Home temp servo for 1.0 Cs MOT	Lv1 above Till driver		
12	Valon	1.5 amplifier stack Lv1 above 1.0 Cs MOT		
13	Till driver v2 for 1.5 Na Raman	Lv1 above 1.5 Na Raman		
14	Home temp servo for 1.5 Na Raman (x2 + adaptor box)	Lv1 above 1.5 Na Raman		
15	New Focus driver for Cs Raman	Lv1 above 1.0 Cs Raman		
16	Benchtop power	Lv1 next to New Focus driver		
17	Shutter driver	Lv1 above 1.0 Cs Raman		
18	New Focus lock box for 1.0 Cs Raman	Lv1 above 1.0 Cs Raman		
19	Timebase drivers for 1.0 Na Raman and Na D2	Lv1 above Thorlabs temp servo		
20	Thorlabs temp servo for 1.0 Na Raman doubler	Lv1 above 1.0 Na seed		
21	Home lock box for Na D2	Lv1 above 1.0 Na seed		
22	MPB Raman laser	Lv1 above 1.0 Na MOT		
23	PLL for Cs Raman	Next to 1.0 Cs Raman		
24	PLL for Cs MOT	Below Cs MOT		

#### 1 Level 2 cloud

• Home power  $(\pm 15,5)$  [1]

#### Turn off check list

- Check On UPS

• Benchtop power 2+1 Chn [2]
Lv2 above 1.0 Cs MOT facing computer
Record voltage/current values:

	Voltage	Current		
Left				
Right				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • Benchtop power 1 Chn [3] Lv2 above 1.0 Cs MOT facing computer Record voltage/current values:

Voltage	Current

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • Benchtop power 1 Chn (x3)

[4]

#### Lv2 above 1.0 Cs MOT facing wall

Record voltage/current values:

	Voltage	Current	UPS	OFF
1				
2				
3				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • Home power (24) (x2) [5] Lv2 above 1.0 Cs Raman

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • Benchtop power for 1.0 Cs Raman [6]

Lv2 above New Focus lock box Record voltage/current values:

Voltage	Current			

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • Benchtop power for 1.0 Na Raman seed [8]

Lv2 above oscilloscope
Record voltage/current values:

Voltage	Current				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

• III High voltage power [7]

Lv2 above timebase drivers Follow home lock box for Na D2 [21].

#### 2 Level 1 clound

• Benchtop power 2+1 Chn for Cs MOT beat lock [9] Lv1 above 1.0 Cs MOT facing computer

Record voltage/current values:

	Voltage	Current		
Left				
Right				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • Till driver v2 for 1.0 Cs MOT [10]

Lv1 above benchtop powere supply Record current values:

Current	

#### Turn off check list

- Turn off
- Unplug power (zip tie to controller)

#### Turn on check list

- Plug in power
- Turn on

### • Home temp servo for 1.0 Cs MOT [11]

Lv1 above Till driver

#### Turn off check list

– Check On UPS



1.5 amplifier stack Lv1 above 1.0 Cs MOT Check with 1.5 to record setting

Frequency	
Amplitude	

#### Turn off check list

– Unplug power from wall

- Plug in power
- Restore settings

•	Off	$\mathbf{T}$	ill	driver	v2	for	1.5	Na	Ra-
	man	[13]							

Lv1 above 1.5 Na Raman

Record current values:

Current

Turn off check list

– Turn off

- Unplug power (zip tie to controller)

Turn on check list

– Plug in power

– Turn on

# • New Focus driver for Cs Raman [15]

Lv1 above 1.0 Cs Raman

Record Temperature, Current, Piezo voltage Temperature is accessible in the menu under system status

F3 Current	
F3 Piezo Voltage	
F3 Temperature	
F4 Current	
F4 Piezo Voltage	
F4 Temperature	

#### Turn off check list

– Turn off output.

– Check plugged into UPS.

#### Turn on check list

- (Turn on)

Home temp servo for 1.5 Na Raman (x2 + adaptor box) [14]
Lv1 above 1.5 Na Raman

Turn off check list

– Check On UPS

### • Benchtop power [16]

Lv1 next to New Focus driver Record voltage/current values:

Voltage	Current

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

• Shutter driver [17]
Lv1 above 1.0 Cs Raman
Switch off/on with the switch on the back.

# • New Focus lock box for 1.0 Cs Raman [18] Lv1 above 1.0 Cs Raman

#### Turn off check list

- Unlock
- Turn down VGA gain
- Power off from the back

#### Turn on check list

- Power on
- Make sure PLL setting is restored
- Turn up VGA gain
- Lock

### • Timebase drivers for 1.0 Na Raman and Na D2 [19]

Lv1 above Thorlabs temp servo Record current and T1 settings

Raman Current	
Raman T1	
D2 Current	
D2 T1	

#### Turn off check list

- Turn off output.
- Check plugged into UPS.

#### Turn on check list

– Turn on

# • Thorlabs temp servo for 1.0 Na Raman doubler [20]

Lv1 above 1.0 Na seed

Record temperature settings

Setpoint	
1	

#### Turn off check list

- Check On UPS

### • Home lock box for Na D2 [21]

Lv1 above 1.0 Na seed

Turn on-off together with the adjustable high voltage power supply above

Record high voltage power supply output

Voltage
---------

#### Turn off check list

- Unlock
- Turn down the voltage of the high voltage power supply to 0
- Unplug high voltage power on lockbox (zip tie)
- Unplug 15V from lockbox (zip tie)
- Turn off high voltage power supply
- Unplug high voltage power supply from wall

#### Turn on check list

- Plug high voltage power supply into the wall
- Turn on high voltage power supply (make sure the output is 0)
- Plug in 15V to lockbox
- Plug in high voltage to lockbox
- Turn up the voltage of the high voltage power supply
- Check HV output

# • MPB Raman laser [22] Lv1 above 1.0 Na MOT

#### Turn off check list

- Power off
- Unplug from wall

#### Turn on check list

- Plug into wall
- Power on

#### 3 Table

## • PLL for Cs Raman [23] Next to 1.0 Cs Raman

#### Turn off check list

- Unplug all signals (zip tie together)
- Unplug 15V and 5V powers (zip tie together)

- Plugin 15V and 5V powers
- Plugin all signals
- Apply settings

#### 4 Floor

# • PLL for Cs MOT [24] Below Cs MOT

#### Turn off check list

- − Unplug all signals (zip tie together)
- Unplug 15V and 5V powers (zip tie together)

- Plugin 15V and 5V powers
- Plugin all signals
- Apply settings

#### Kenneth

For shutting off, go through the checklist before turning off the computers.

- Turn off equipments
- Turn off computers

Wait before everyone else finish their equipment shutoff before turning off the computers. For powering up, first turn the computers on before turning on the equipments.

- Turn on computers
- Turn on equipments
- Check computer connection to devices

Also wait before everyone else finish their corresponding steps before starting the next one.

#### Record settings before turning things off.

All the lasers should be turned off.

All the home built  $\pm 15$  V, 5V power supplies should be on UPS. 24 V ones are mainly powering amplifiers and can be shut down.

All locks should be disabled.

All temperature controllers should be on UPS and left on.

For each equipment, check if it is in used. Then check the power and make sure if is/should be on UPS. Make sure the "On UPS" label is correct.

## Equipments are mostly located around the laser table facing the machine table.

### Equipment list

ID	Equipment	Location	
	Laser table side facing machine table top $\rightarrow$ down, north side $\rightarrow$ east side		
1	Benchtop power 1 Chn	Lv2 northwest corner	
2	Home power (24)	Lv2 above 623 nm box	
3	Benchtop power 2+1 Chn	Lv2 above D1 doubler beam path	
4	Benchtop power 1 Chn	Lv2 left of 2 Chn benchtop power	
5	FPGA box for 1.0	Lv2	
6	Software radio	Lv2	
7	High voltage power	Lv2 next to software radio	

8	Benchtop power 2+1 Chn	Lv2 next to radio computer
9	Radio computer	Lv2 northeast corner
10	Home temp servo	Lv1 next to oscilloscope
11	Timebase 1038 control box	Lv1 above 623 nm box
12	Till driver v1 (unused)	Lv1 above D1 double beam path
13	Thorlabs temp servo for Na D1 seed	Lv1 above Na D1
14	Thorlabs current for Na D1 seed	Lv1 above Na D1
15	Home temp servo	Lv1 above 976 nm ECDL (Mango)
16	Thorlabs temp servo for Na D1 doubler	Lv1 above Na D1
17	Home lock box for Na D1	Lv1 above Na saturated absorption beam path
18	Shutter driver	Lv1 next to variac
19	Variac for Na cell	Lv1 above 1.0 Na MOT
20	NKT amplifier	Floor northwest corner
21	SFG temp servo	Floor on NKT amplifier
22	Home temp servo	Floor on NKT amplifier
23	Thorlabs current	Floor on NKT amplifier
24	ALS fiber laser	Floor middle of table
25	671 EOM stack	Floor northeast corner

#### 1 Level 2 cloud

### • Benchtop power 1 Chn [1]

#### Lv2 northwest corner

Record voltage/current values:

Voltage	Current

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • Home power (24) [2] Lv2 above 623 nm box

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • Benchtop power 2+1 Chn [3]

#### Lv2 above D1 doubler beam path

Record voltage/current values:

	Voltage	Current
Left		
Right		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • Benchtop power 1 Chn [4]

Lv2 left of 2 Chn benchtop power Record voltage/current values:

Voltage	Current

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • FPGA box for 1.0 [5]

#### Lv2

#### Turn off check list

- (Do following three steps quickly)
- Turn off the box by flipping the switch on the front side of the box
- Unplug the high voltage (48V) power supply (circular plug) on the front side of the box. Then unplug the power supply itself fromm the outlet.
- Unplug the 12V board power supply. (L.T.E. one with a rectangular connector in the front)
- Unplug USB power on the front side near the fan.
- Zip tie the three unpluged powers to the computer control box.

#### Turn on check list

- Connect to 3.5G Windfrek clock generator from Tweezer1. Set frequency to 3.5G and amplitude to max.
- (Remove power connection zip ties and) Plug in USB power supply.
- (Do following three steps quickly)
- Plug in the 12V board power supply.
- Plug in high voltage power supply.
- Turn on the power switch of the box. The fans should start spinning.

### • Software radio [6]

Turn off/on with the switch on the front.

# • High voltage power [7] Lv2 next to software radio

Follow home lock box for Na D1 [17].

### • Benchtop power 2+1 Chn [8]

#### Lv2 next to radio computer

Record voltage/current values:

	Voltage	Current
Left		
Right		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • Radio computer [9]

Ask Yichao

#### 2 Level 1 cloud

• Home temp servo [10]

#### Turn off check list

– Check On UPS

Timebase 1038 control box  [11] Lv1 above 623 nm box  Should be on UPS and Off.  Turn off check list   Turn on check list	Thorlabs current for Na D1 seed [14] Lv1 above Na D1 Record current settings  Current  Turn off check list  - Turn off - Unplug power  Turn on check list  - Plug in power  - Check current setting  - Turn on
Till driver v1 (unused) [12] Lv1 above D1 double beam path Turn off check list  Check On UPS	Home temp servo [15] Lv1 above 976 nm ECDL (Mango)  Turn off check list  Check On UPS
Thorlabs temp servo for Na D1 seed [13] Lv1 above Na D1 Record temperature settings  Setpoint  Turn off check list  - Check On UPS	Thorlabs temp servo for Na D1 doubler [16] Lv1 above Na D1 Record temperature settings  Setpoint  Turn off check list  - Check On UPS

# • Home lock box for Na D1

Lv1 above Na saturated absorption beam path Turn on-off together with the fixed high voltage power supply above

#### Turn off check list

- Unlock
- Unplug high voltage power on lockbox (zip tie)
- Unplug 15V from lockbox (zip tie)
- Unplug high voltage power supply from walls

#### Turn on check list

- Plug high voltage power supply into the wall
- Plug in 15V to lockbox
- Plug in high voltage to lockbox
- Check HV output

### • Variac for Na cell [19]

Lv1 above 1.0 Na MOT

Record temperature and voltage settings

Voltage	
Temperature 1	
Temperature 2	

Check if we can use the backup power. (Through the UPS)  $\bf Turn\ off\ check\ list$ 

- Turn voltage to 0
- Switch off
- Unplug from wall

#### Turn on check list

- Plug into wall
- Switch on
- Turn the voltage up **SLOWLY** and monitor the temperature at the same time.

# • Shutter driver [18] Lv1 next to variac

Switch off/on with the switch on the back.

0	Floor
.3	HIOOT
v	1 1001

### • NKT amplifier [20]

Floor northwest corner

Record settings,

ACC	
Pin	

#### Turn off check list

- Power to 0
- Turn off
- Key off
- Unplug power

#### Turn on check list

- Plug in power
- Key on
- Turn on
- Restore settings

### • SFG temp servo [21]

Floor on NKT amplifier Record temperature settings

ı		
ı		
ı	G	
ı	Setpoint	
ı	Desponie	

#### Turn off check list

– Check On UPS

# • Home temp servo [22] Floor on NKT amplifier

#### Turn off check list

- Check On UPS

### Thorlabs current [23]

Floor on NKT amplifier Record current settings

Current

#### Turn off check list

- Turn off
- Unplug power

- Plug in power
- Check current setting
- Turn on

# ALS fiber laser [24] Floor middle of table

Record settings,

Power	
Current	

#### Turn off check list

- Power to 0
- Turn off
- Key off
- Unplug power

#### Turn on check list

- Plug in power
- Key on
- Turn on
- Restore settings



Floor northeast corner

Disassemble?

#### Lewis

For shutting off, go through the checklist before turning off the computers.

- Turn off equipments
- Turn off computers

Wait before everyone else finish their equipment shutoff before turning off the computers. For powering up, first turn the computers on before turning on the equipments.

- Turn on computers
- Turn on equipments
- Check computer connection to devices

Also wait before everyone else finish their corresponding steps before starting the next one.

#### Record settings before turning things off.

All the lasers should be turned off.

All the home built  $\pm 15$  V, 5V power supplies should be on UPS. 24 V ones are mainly powering amplifiers and can be shut down.

All locks should be disabled.

All temperature controllers should be on UPS and left on.

For each equipment, check if it is in used. Then check the power and make sure if is/should be on UPS. Make sure the "On UPS" label is correct.

## Equipments are mostly located around the machine table facing the laser table.

### Equipment list

ID	Equipment	Location			
	Machine table facing the laser table top $\rightarrow$ down, west side $\rightarrow$ south side				
1	Keithley power (x2)	Lv2 northwest corner			
2	Benchtop power 2+1 Chn	Lv2 facing computer			
3	Benchtop power 2 Chn	Lv2 above 1.5 Cs MOT/RP current			
4	Home power (24) x2	Lv2 next to benchtop power			
5	Benchtop power	Lv2 on two home power (24)			
6	Home power $(\pm 15,5)$	Lv2 middle			
7	Benchtop power (x2, 1 Chn & 2+1 Chn) for uWave Amp	Lv2 next to home power $(\pm 15,5)$			

8	Benchtop power 1 Chn	Lv2 above uWave breadboard
9	Benchtop power (x2, 1 Chn & 2+1 Chn)	Lv2 southeast corner
10	Function generator for Na switching	Lv1 facing computer
11	Till driver v1 for $1.5~\mathrm{Cs}~\mathrm{RP}$	Lv1 above 1.5 Cs
12	Till driver v2 for $1.5~\mathrm{Cs}~\mathrm{MOT}$	Lv1 above 1.5 Cs
13	Home temp servo (x2) for 1.5 Cs MOT/RP	Lv1 above 1.5 Cs
14	Home lock box for 1.5 Cs	Lv1 above 1.5 Cs
15	Novatec for 1.5 Cs lock	Lv1 above 1.5 Cs
16	Benchtop power 1 Chn (x2)	Lv1 above 1.5 Na
17	High voltage power	Lv1 next to 2 benchtop power
18	Function generator for 1.0 MOT piezo	Lv1
19	Temperature and pump controllers for cavity	Lv1 next to function generator
20	Power (for Na uWave)	Lv1 next to uWave breadboard
21	SAES pump	Lv1 above 1.0 chamber
22	PLL for 1.5 Cs MOT	Table near 1.5 Cs MOT
23	Controllers for TiSapph (Mercury)	Floor northwest corner
24	Controllers for TiSapph (Venus)	Floor southwest corner
25	Water chiller for uWave breadboard	Floor below 1.5 Na MOT beam path

#### 1 Level 2 cloud

### • Keithley power (x2) [1]

#### Lv2 northwest corner

Record voltage/current values:

	Voltage	Current	UPS	OFF
1				
2				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### Benchtop power 2+1 Chn [2]

#### Lv2 facing computer

Record voltage/current values:

	Voltage	Current
Left		
Right		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • Benchtop power 2 Chn [3]

Lv2 above 1.5 Cs MOT/RP current Record voltage/current values:

	Voltage	Current
Left		
Right		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • Home power (24) x2 [4] Lv2 next to benchtop power

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • Benchtop power [5]

Lv2 on two home power (24) Record voltage/current values:

Voltage	Current

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- $\square$  (if turn off) Plug the load back in
- Home power  $(\pm 15,5)$  [6]

#### Turn off check list

- Check On UPS
- Benchtop power (x2, 1 Chn & 2+1 Chn) for uWave Amp [7]
  Lv2 next to home power (±15,5)

Off

### • Benchtop power 1 Chn [8]

Lv2 above uWave breadboard Record voltage/current values:

Voltage	Current		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • Benchtop power (x2, 1 Chn & 2+1 Chn) [9]

Lv2 southeast corner

Record voltage/current values:

	Voltage	Current	UPS	OFF
1				
2 L				
2 R				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

#### 2 Level 1 cloud

# • Function generator for Na switching [10]

Lv1 facing computer

Record High, Low, Frequency, Phase, Duty cycle

No. 1 High V	
No. 1 Low V	
No. 1 Frequency	
No. 1 Phase	
No. 1 Duty cycle	
No. 2 High V	
No. 2 Low V	
No. 2 Frequency	
No. 2 Phase	
No. 2 Duty cycle	

#### Turn off check list

- Check On UPS

# • Till driver v1 for 1.5 Cs RP

Lv1 above 1.5 Cs

Record current values:

Current

Turn off check list

- Turn off
- Unplug power (zip tie to controller)

Turn on check list

- Plug in power
- Turn on
- Till driver v2 for 1.5 Cs MOT [12]

Lv1 above 1.5 Cs

Record current values:

Current

Turn off check list

- Turn off
- Unplug power (zip tie to controller)

Turn on check list

- Plug in power
- Turn on
- Home temp servo (x2) for 1.5 Cs MOT/RP [13]

Lv1 above 1.5 Cs

Turn off check list

- Check On UPS

• 1.5 Home lock box for 1.5 Cs

Lv1 above 1.5 Cs

Turn off check list

- Unlock
- Unplug 15V from lockbox (zip tie)

Turn on check list

– Plug in 15V to lockbox

• Novatec for 1.5 Cs lock [15]
Lv1 above 1.5 Cs

Turn off check list

- Save settings

Turn on check list

- Restore settings

# • Benchtop power 1 Chn (x2) [16]

Lv1 above 1.5 Na

Record voltage/current values:

	Voltage	Current	UPS	OFF
1				
2				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • High voltage power [17] Lv1 next to 2 benchtop power Follow function generator for 1.0 MOT piezo [18].

# • Function generator for 1.0 MOT piezo [18]

Turn on-off together with the HV amplifiers for 1.0 piezo mirrors above and the high voltage power supply on the left

Turn off check list

- Turn off the function generator
- Unplug function generator output (zip tie and label)
- Unplug HV power supply from the wall

#### Turn on check list

- Plugin HV power supply to the wall
- Plugin function generator output
- Turn on the function generator

# • Temperature and pump controllers for cavity [19] Lv1 next to function generator

#### Turn off check list

– Check On UPS

## • Power (for Na uWave) [20] Lv1 next to uWave breadboard

Unused.

# • SAES pump [21] Lv1 above 1.0 chamber

#### Turn off check list

– Check On UPS

#### 3 Table

• PLL for 1.5 Cs MOT [22]
Table near 1.5 Cs MOT

#### Turn off check list

- Unplug all signals (zip tie together)
- Unplug 15V and 5V powers (zip tie together)

#### Turn on check list

- Plugin 15V and 5V powers
- Plugin all signals
- Apply settings

#### 4 Floor

• Controllers for TiSapph (Mercury) [23]
Floor northwest corner

Floor northwest corner Record settings,

Pump power	
Chiller temperature	

#### Turn off check list

- Standby pump
- Pump off, shutter close, unplug from wall
- Turn off SolsTis controller (Lv1 facing KRb), unplug from wall
- Turn off Chiller, unplug from wall

#### Turn on check list

- Turn on chiller, restore setting.
- Turn on SolsTis controller.
- Check wavelength setting (700 nm)
- Pump controller on, pump standby.
- Pump shutter open. Pump set to original power.

# • Controllers for TiSapph (Venus) [24]

Floor southwest corner Record settings,

Chiller temperature	
Pump power setpoint	
Pump power reading	
Wavelength setting	

#### Turn off check list

- Pump power set to 0, shutter close.
- Pump cool down (takes a long time)
- Pump off, unplug from wall
- Turn off SolsTis controller, unplug from wall
- Turn off Chiller, unplug from wall

#### Turn on check list

- Turn on chiller, restore setting.
- Turn on SolsTis controller.
- Check wavelength setting
- Pump controller on, pump warm up (takes a long time).
- Pump set to original power.

# • Water chiller for uWave breadboard [25] Floor below 1.5 Na MOT beam path

Turn off.

#### Jessie

For shutting off, go through the checklist before turning off the computers.

- Turn off equipments
- Turn off computers

Wait before everyone else finish their equipment shutoff before turning off the computers. For powering up, first turn the computers on before turning on the equipments.

- Turn on computers
- Turn on equipments
- Check computer connection to devices

Also wait before everyone else finish their corresponding steps before starting the next one.

#### Record settings before turning things off.

All the lasers should be turned off.

All the home built  $\pm 15$  V, 5V power supplies should be on UPS. 24 V ones are mainly powering amplifiers and can be shut down.

All locks should be disabled.

All temperature controllers should be on UPS and left on.

For each equipment, check if it is in used. Then check the power and make sure if is/should be on UPS. Make sure the "On UPS" label is correct.

#### Equipments are mostly located around machine table facing KRb.

### Equipment list

ID	Equipment	Location		
	Machine table facing KRb top $\rightarrow$ down, north side			
1	Highest (voltage) power in the lab	Zip tied to you know where it is		
2	Benchtop power 1 Chn	Lv2 northwest corner facing KRb		
3	Home power $(\pm 15,5)$	Lv2 above 1.5 Cs Raman		
4	Benchtop power 1 Chn (x2)	Lv2 on home power $(\pm 15,5)$		
5	High voltage power	Lv2 between home power ( $\pm 15,5$ ) and FPGA box		
6	FPGA box for 1.5	Lv2		
7	Valon	Lv2 left of FPGA box		

8	Home temp servo x3 for STIRAP	Lv2 left of FPGA box
9	Benchtop power (x9, 2+1 Chn x2 & 2 Chn x1 & 1 Chn x6)	Lv2
10	Home power (24)	Lv2 above 1.0 tweezer breadboard
11	Home power $(\pm 15,5)$	Lv2 on home power (24)
12	Benchtop power 2 Chn	Lv2 northeast corner
13	Keithley multimeter (x2)	Lv2 northeast corner
14	Shutter driver	Lv1 on M2 controller
15	Function generator	Lv1 on shutter driver
16	Home temp servo	Lv1 left of M2 controller
17	Home servo x2 for 1.5 tweezers	Lv1 left of home temp servo
18	Till driver v2	Lv1 above 1.5 Cs Raman beam path
19	Delay generator for 1.5 Na	Lv1 below 1.5 output panel
20	Rack for STIRAP	Lv1 above STIRAP
21	Till driver (x2, v1 & v2) for STIRAP lasers	Lv1 above rack for 1.5
22	Vincent servo x2 for STIRAP lasers	Lv1
23	Home servo x2	Lv1
24	Current servo for compensation coil	Lv1 below home power (24)
25	Greiner servo	Lv1 northeast corner
26	Home temp servo x2 for 1.5 Cs Raman	Floor rack Lv2
27	Piezo driver	Floor rack Lv2 on home temp servo x2
28	Benchtop power 1 Chn	Floor rack Lv2 left of uWave foam
29	Till driver v2 x2 for 1.5 Cs Raman	Floor rack Lv2
30	Home servo for 1.5 Cs Raman phase lock	Floor rack Lv2 above Till drivers
31	Benchtop power 2+1 Chn	Floor rack Lv1
32	Benchtop power 1 Chn	Floor rack Lv1

33	Chiller for 1.0 camera	Floor
34	Benchtop power 1 Chn x 3	Floor STIRAP rack Lv2

#### 1 Level 3 cloud

• Highest (voltage) power in the lab [1]
Zip tied to you know where it is...

#### 2 Level 2 cloud

• Benchtop power 1 Chn [2]
Lv2 northwest corner facing KRb
Record voltage/current values:

Voltage	Current

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- − (if turn off) Plug the load back in

### • Benchtop power 1 Chn (x2)

[4]

2

 $\overline{\text{Lv2}}$  on home power (±15,5) Record voltage/current values:

	Voltage	Current	UPS	OFF
1				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- ☐ (if turn off) Plug the load back in

# • Home power $(\pm 15,5)$ [3]

#### Turn off check list

– Check On UPS

# • High voltage power [5] Lv2 between home power (±15,5) and FPGA box

Check connection and power off.

### • FPGA box for 1.5 [6]

#### Lv2

Save 3.5 GHz amplitude

Amplitude
-----------

#### Turn off check list

- Save 3.5 GHz amplitude.
- (Do following three steps quickly)
- Turn off the box by flipping the switch on the front side of the box
- Unplug the high voltage (48V) power supply (circular plug) on the left side of the box. Then unplug the power supply itself fromm the outlet.
- Unplug the 12V board power supply. (L.T.E. one with a rectangular connector in the front)
- Unplug USB power on the front side near the fan.
- Zip tie the three unpluged powers to the computer control box.

#### Turn on check list

- Set 3.5 GHz frequency and amplitude.
- (Remove power connection zip ties and) Plug in USB power supply.
- (Do following three steps quickly)
- Plug in the 12V board power supply.
- Plug in high voltage power supply. The fans should start spinning.
- Turn on the power switch of the box.

# • Valon [7] Lv2 left of FPGA box

Save frequency and amplitude

	Frequency	Amplitude
Ch 1		
Ch 2		

#### Turn off check list

– Unplug power from wall

#### Turn on check list

- Plug in power
- Restore settings

Home temp servo x3 for STI-RAP [8] Lv2 left of FPGA box

#### Turn off check list

- Check On UPS

# Benchtop power (x9, 2+1 Chn x2 & 2 Chn x1 & 1 Chn x6) $_{\rm Lv2}^{\rm [9]}$

Record voltage/current values:

	Voltage	Current	UPS	OFF	
1 Channel x 6					
1					
2					
3					
4					
5					
6					
	2 C	hannel x 1			
7 L					
7 R					
2+1 Channel x 2					
8 L					
8 R					
9 L					
9 R					

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged 6
- (if turn off) Plug the load back in

### Home power (24) [10] Lv2 above 1.0 tweezer breadboard

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged
- (if turn off) Plug the load back in

- Home power  $(\pm 15,5)$  [11] Lv2 on home power (24) Turn off check list
  - Check On UPS

### • Benchtop power 2 Chn [12]

Lv2 northeast corner

Record voltage/current values:

	Voltage	Current
Left		
Right		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

Keithley multimeter (x2) [13]

Lv2 northeast corner

Ask Kenneth

#### 3 Level 1 cloud

• Shutter driver [14]

Lv1 on M2 controller

Switch off/on with the switch on the back.

# • Function generator [15] Lv1 on shutter driver

Unused. Remove.

- Home temp servo [16]

  Lv1 left of M2 controller

  Turn off check list
  - Check On UPS

# Home servo x2 for 1.5 tweezers [17] Lv1 left of home temp servo

#### Turn off check list


- Unplug power
- Tie power to box

#### Turn on check list

– Plug power

## Till driver v2 [18]

Lv1 above 1.5 Cs Raman beam path Record current values:

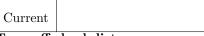

#### Turn off check list

- Turn off
- Unplug power (zip tie to controller)

- Plug in power
- Turn on

Delay generator for 1.5 Na  [19] Lv1 below 1.5 output panel Record current values:	Till driver (x2, v1 & v2) for STIRAP lasers [21] Lv1 above rack for 1.5 Record current values:
	Current 1  Current 2  Turn off check list  - Turn off  - Unplug power (zip tie to controller)
Turn off check list  -	Turn on check list  - Plug in power  - Turn on
Turn on check list  -	
• Rack for STIRAP [20] Lv1 above STIRAP Ask Lewis	$\begin{array}{c ccccc} \bullet & & & & & & & & & & & & & & & & & & $
Turn off check list  -	Turn off check list  -
Turn on check list	Turn on check list

# • Home servo x2 [23]

Turn off check list

- Unplug power
- Tie power to box

Turn on check list

- Plug power
- Current servo for compensation coil [24]

Lv1 below home power (24) Ask Lewis

Turn off check list



Turn on check list



• Greiner servo [25]

Turn off check list

- Unplug power
- Tie power to box

Turn on check list

– Plug power

- 4 Floor
- Home temp servo x2 for 1.5 Cs Raman [26]

Turn off check list

- Check On UPS

• Piezo driver [27]
Floor rack Lv2 on home temp servo x2

Check high voltage off/on.

• Benchtop power 1 Chn [28]
Floor rack Lv2 left of uWave foam
Record voltage/current values:

Voltage	Current	

Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply **Do NOT do this with loads pluged** in
- (if turn off) Plug the load back in

### • Till driver v2 x2 for 1.5 Cs Raman [29]

Floor rack Lv2

Record current values:

Current 1	
Current 2	

#### Turn off check list

- Turn off
- Unplug power (zip tie to controller)

#### Turn on check list

- Plug in power
- Turn on

# • Benchtop power 2+1 Chm

Floor rack Lv1

Record voltage/current values:

	Voltage	Current
Left		
Right		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

# • Home servo for 1.5 Cs Raman phase lock [30] Floor rack Lv2 above Till drivers

#### Turn off check list

- Unplug power
- Tie power to box

#### Turn on check list

– Plug power

### • Benchtop power 1 Chn [32]

#### Floor rack Lv1

Record voltage/current values:

Voltage	Current

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- − (if turn off) Plug the load back in

### • Chiller for 1.0 camera [33]

#### Floor

Record temperature setting,

Setpoint	

#### Turn off check list

- Turn off camera
- Turn off chiller
- Unplug power

#### Turn on check list

- Plug power
- Turn on chiller
- Turn on camera

# • Benchtop power 1 Chn x 3 [34]

#### Floor STIRAP rack Lv2

Record voltage/current values:

	Voltage	Current	UPS	OFF
1				
2				
3				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

#### Will

For shutting off, go through the checklist before turning off the computers.

- Turn off equipments
- Turn off computers

Wait before everyone else finish their equipment shutoff before turning off the computers. For powering up, first turn the computers on before turning on the equipments.

- Turn on computers
- Turn on equipments
- Check computer connection to devices

Also wait before everyone else finish their corresponding steps before starting the next one.

#### Record settings before turning things off.

All the lasers should be turned off.

All the home built  $\pm 15$  V, 5V power supplies should be on UPS. 24 V ones are mainly powering amplifiers and can be shut down.

All locks should be disabled.

All temperature controllers should be on UPS and left on.

For each equipment, check if it is in used. Then check the power and make sure if is/should be on UPS. Make sure the "On UPS" label is correct.

# Equipments are mostly located around 1.5 table. The power supply on the second floor may need to be shut down first.

### Equipment list

ID	Equipment	Location	
	1.5 Rack		
1	Home power (24)	Bottom level	
2	Benchtop power 2 Chn	Above home power (24)	
3	SAES pump controller	Middle left	
4	Benchtop power 1 Chn x2 for dispensers	Middle next to pump controller	
5	Benchtop power 2 Chn for 1.5 coil	Above oscilloscopes	
6	Home power $(\pm 15,5)$	Top level	

7	Function generator for 1.5 MOT piezo	Top level on high voltage Amp		
	Below 1.5 table			
8	IGBT stack	Floor		
9	Chiller for 1.5 camera	Floor behind IGBT		
10 Benchtop power 1 Chn x 4		Floor between table and drawers		
	Other			
11	Feshbach power	HVAC room		
12	Cooling water			

#### 1 1.5 Rack

• Image: Home power (24) [1]

Bottom level Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- ☐ (if turn off) Plug the load back in

### • Enchtop power 2 Chn [2]

Above home power (24)

Record voltage/current values:

	Voltage	Current
Left		
Right		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- − (if turn off) Plug the load back in

# • SAES pump controller [3]

#### Turn off check list

– Check On UPS

# • Benchtop power 1 Chn x2 for dispensers [4] Middle next to pump controller

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- ☐ (if turn off) Plug the load back in

# • Benchtop power 2 Chn for 1.5 coil [5]

Above oscilloscopes

Record voltage/current values:

	Voltage	Current
Left		
Right		

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

• Home power  $(\pm 15,5)$  [6]

Turn off check list

– Check On UPS

#### Function generator for 1.5 Chiller for 1.5 camera [9] Floor behind IGBT MOT piezo [7] Floor Top level on high voltage Amp Record temperature setting, Floor Record settings, Setpoint Frequency Turn off check list Amplitude – Turn off camera Turn off chiller Turn off check list Unplug power - Turn off – Unplug power Turn on check list Plug power Turn on check list Turn on chiller Plug power

Turn on camera

#### 2 Below 1.5 table

Turn on

• IGBT stack [8]

#### Turn off check list

- Make sure Feshbach current is off
- Turn off/unplug Feshbach servo box

- Plug in/turn on Feshbach servo box
- Make sure IGBTs are still working

### Benchtop power 1 Chn x 4

[10]

Floor between table and drawers

Record voltage/current values:

	Voltage	Current	UPS	OFF
1				
2				
3				
4				

#### Turn off check list

- Decide (On UPS / Turn off)
- (if on UPS) Check On UPS
- (if turn off) Unplug load
- (if turn off) Label and zip tie loads to the power supply
- (if turn off) Turn off and unplug from wall

#### Turn on check list

- (if turn off) Plug into the wall and turn on the power supply Do NOT do this with loads pluged in
- (if turn off) Plug the load back in

### • Feshbach power [11]

11 1/10 100111		
	Voltage	
	Current	
Record settings,		

#### Turn off check list

- Switch off.
- L

#### Turn on check list

- \_ \_
- \_
- Make sure water cooling is on.
- Switch on.
- Restore settings.

### • Cooling water [12]

- Turn off check list
  - Make sure coil power and IGBTs are off.
  - Switch off.

#### Turn on check list

– Switch on before turning on coil power.